




THE POSTURAL LINE

The Walker-View is a fully equipped treadmill with an embedded PC platform capable of simultaneously assessing both gait analysis and the movement of all parts of the body. The 47 inch monitor, positioned at the front of the Walker-View system, is a digital mirror providing the client with powerful feedback on controlling and analyzing their posture when in movement.

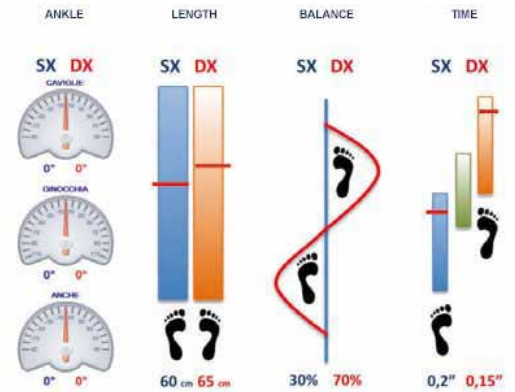

WALKER-VIEW




- Standard Feature
- Optional Feature

USER INTERFACE

SPECIFICATIONS

	WALKER-VIEW
47" LCD Feedback Display	●
15" PC Control Touchscreen	●
Processor Frequency	2.5 GHz
RAM Memory	4GB
Wireless Keyboard & Mouse	●
Printable Reports (printer not included)	●
TecnoBody Management System	●
Running Surface Dimensions	20" W x 63" L (51 x 160 cm)
User Weight Capacity	330 lb.
Unit Weight	617 lb.
Width	35.4" (90 cm)
Length	108.3" (275 cm)
Height	80.7" (205 cm)
0-15% Incline	●
0-12.5 MPH	●
Reverse 0-5 MPH	●
3D Camera & Force Platform	●
208/230 Vac 20 Amp	●
Warranty	1 Year All Components 1 Year Labor Lifetime Software

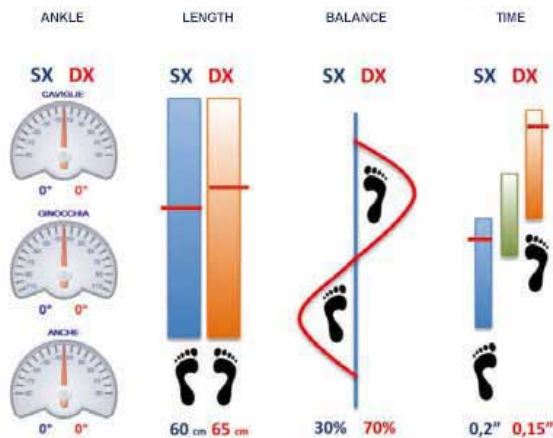

SOFTWARE

- Test modules for Gait Analysis, Posture Control, and Aerobic Endurance
- Integrated 3D camera to analyze movement and posture
- Virtual reality games for functional and interactive walking & running



FUNCTIONALITY

- Embedded PC platform monitors and records step length, step speed, step symmetry, and range of motion of the hips, knees, and ankles
- The Walker-View treadmill is equipped with a sensitive load cell belt allowing the assessment of walking/running stance and corrects dynamic parameters in real time
- Ability to connect to TecnoBody Management System by TecnoBody Key



FOR ADDITIONAL INFORMATION ABOUT THIS PRODUCT VISIT WWW.WOODWAY.COM OR CALL 800.WOODWAY TO TALK TO A REPRESENTATIVE.